

Electronic Pulse Massager

MODEL HV-F128

With 12 Auto Therapy Modes and 5 Special Modes to choose from, the HV-F128 Electronic Pulse Massager provides the perfect level of relief for muscular ache.

LONG-LIFE PADS

- Durable
- Washable

SOOTHE & RELIEVE THERAPY PROGRAM

Helps relieve muscular ache and stiffness

12 | Auto Therapy Modes
5 | Special Modes



All for Healthcare

OMRON

Electronic Pulse Massager

MODEL HV-F128

Enjoy complete therapy effect

- > Relieves muscular aches and numbness
- > Accelerates recovery from fatigue
- > Improves blood circulation
- > Relaxes stiff shoulder and muscle spasms

12 AUTO THERAPY MODES

Mode Selections of Soothe and Relieve for specific regions with massage program selection.

SOOTHE:

Therapy program for muscle stiffness. Low frequency stimulation for effective relief of stiff muscles and improves blood circulation.

RELIEVE:

Therapy program for muscle aches. High frequency stimulation at 1200Hz for effective relief of muscle aches at the joints and lower back.

6 REGION PROGRAM:

Automatic therapy program for specific areas . SHOULDERS . ARMS . JOINTS . WAIST . LEGS . SOLES.

4 MASSAGE PROGRAM:

Choose from any of these 4 programs . TAP . KNEAD . PUSH . RUB.

5 SPECIAL MODES

Use this mode to adjust the therapy intensity required. Combine with Auto Therapy Modes for better effect.

HYPER MODE:

Gradual increase of intensity. Effective for quick relief of muscular ache and spasms.

SOFT MODE:

Gentle stimulation. Effective for recovery from fatigue and muscle aches.

REPEAT MODE:




Repeats your favourite program mode.

POINT & WIDE:

For shoulder or small area- select POINT. For lower back or larger area- select WIDE.

'R' & 'W' MODE:

Adjust therapy intensity on left and right pads as desired.

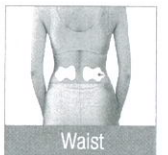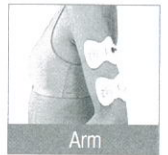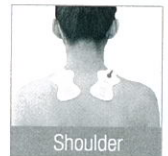

Specifications:

Model	HV-F128
Power Supply	2 x AAA alkaline batteries
Battery Life	Approx. 3 months, when used once a day for 15 minutes Approx. 2 months, when used once a day at HYPER mode for 15 minutes
Frequency Generated	Approx. 1 – 1200Hz
Power Consumption	Approx. 85mA
Unit Dimension	60 (W) x 199 (H) x 21 (D) mm
Unit Weight	Approx. 155g (including batteries)

Warning! This Pulse Massager HV-F128 must not be used in combination with the following medical devices:

- 1) Any internally transplanted electronic devices e.g. pacemakers
- 2) Any electronic life support equipment such as respiratory ventilators
- 3) Any other electronic medical devices attached to the body e.g. ECG machines

This may cause malfunctioning of these devices.

OMRON HEALTHCARE Co.,Ltd
53, Kunotsubo, Terado-cho,
Muko, Kyoto 617-0002 Japan

Asia Pacific HQ
OMRON HEALTHCARE SINGAPORE PTE LTD
438A Alexandra Road, #05-05/08,
Alexandra Technopark, Singapore 119967
www.omron-healthcare.com.sg

Distributed by:

EC REP

OMRON HEALTHCARE EUROPE B.V
Scorpius 33, 2132 LR Hoofddorp
THE NETHERLANDS

CE0197